


# NGS-PRONTO® HIGH-RESOLUTION NANOPORE HLA TYPING

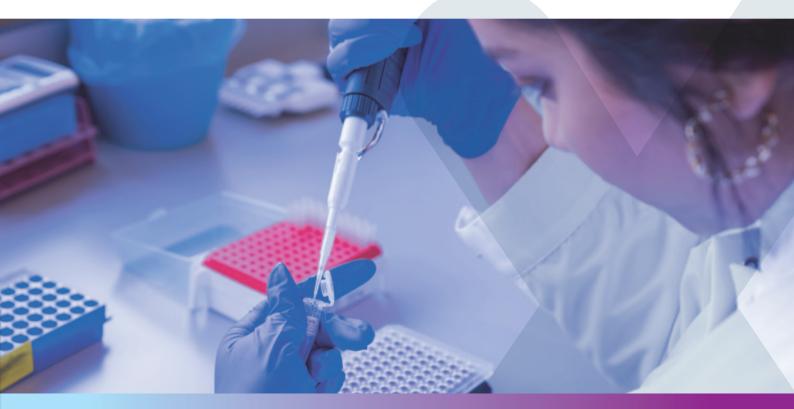

experts in transplant diagnostics

#### Oxford Nanopore Sequencing

11 loci whole gene4 to 96 samples per runAmplification & library preparation

## **THE FUTURE OF NANOPORE HLA TYPING** 12 samples DNA to results in less than 8 hours




Discover our next-level innovation in HLA typing, specifically designed for use with Oxford Nanopore Technologies sequencing devices.

NGS-Pronto offers:

- Easy-to-use amplification and library preparation reagents
- 96 unique barcodes
- An efficient workflow designed for speed and accuracy

NGS-Pronto delivers precise and reliable results, setting a new standard in high-resolution HLA typing.

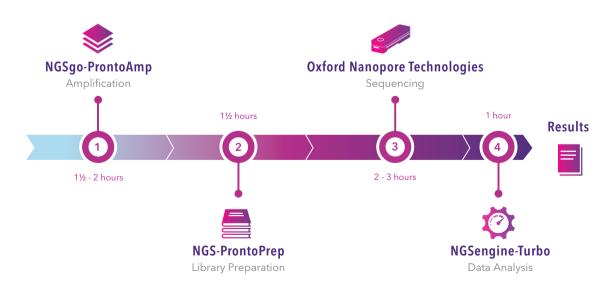

Turnaround time from DNA to results for 12 samples

### AMPLIFICATION & LIBRARY PREPARATION MADE SIMPLE

#### Amplification


NGSgo-ProntoAmp offers amplification of 11 HLA loci using a whole-gene approach, resulting in at least a 3 field resolution.

#### **Library Preparation**

NGS-ProntoPrep offers a complete barcoded library preparation for sequencing on Oxford Nanopore Technologies devices. The procedure only takes 1.5-2 hours and is suitable for 4 to 96 samples. Since there is no fragmentation, long intact amplicons are generated that enable full phasing.

Library preparation is performed at room temperature, allowing for a flexible workflow.



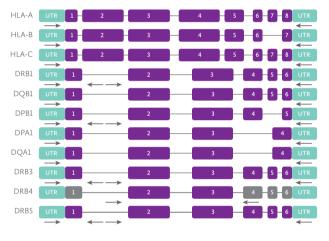

NGSgo®-ProntoAmp offers a whole gene approach

NGS-Pronto® complete kit

### FAST & AFFORDABLE NANOPORE SEQUENCING


NGS-Pronto is compatible with R10.4.1 flow cells and can be sequenced on MinION Mk1B and GridION sequencers from Oxford Nanopore Technologies, delivering rapid results in as little as 15 minutes per sample.



MinION Flow Cell

# DATA ANALYSIS IN MINUTES

Streamline your analysis with NGSengine-Turbo, our user-friendly software tailored for Oxford Nanopore Technologies data. With a fast analysis time of approximately 5 minutes per sample, you'll quickly have the results you need.

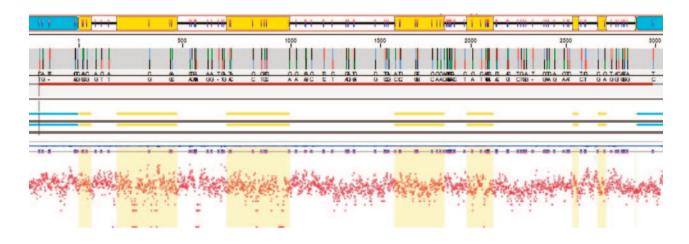

Representative data analyzed by NGSengine-Turbo

# CHALLENGE US WITH YOUR QUESTIONS

Our Support team can assist you with product questions, practical matters, data interpretation, and more.

support@gendx.com +31 30 252 37 99 +1 312 815 5006 (USA)

#### experts in transplant diagnostics



Alexander Numan Building Yalelaan 48 3584 CM Utrecht the Netherlands Phone: +31 30 2523799 Cumberland Metro Office Park 5521 N. Cumberland Avenue Suite 1116 Chicago, IL 60656 USA Phone: +1 312 815 5006 info@gendx.com www.GenDx.com GenDx is a registered trade name of Genome Diagnostics B.V. All other trade names are property of their respective owners, more info www.GenDx.com. M24-001, V1, 2024-08